

Reflex™



Reflex™



- Non-adjustable cushion with fixed amount of immersion
- Simple, no maintenance
- For clients who are at moderate risk for skin breakdown with mild positioning needs
- HCPCS Code E2603 Skin Protection Wheelchair Seat Cushion, non-adjustable

Reflex™ Foam


Two stiffnesses of foam:



- Medium stiffness foam creates stable base and area for loading thighs
- Firmest stiffness foam provides LE positioning as well as posterior and lateral pelvic stability

Reflex™ Foam

Two stiffnesses of foam:

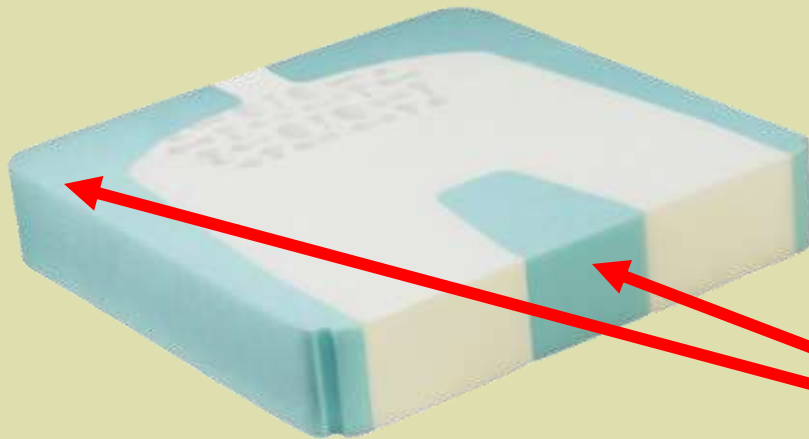


- Medium stiffness foam creates stable base and area for loading thighs
- Firmest stiffness foam provides LE positioning as well as posterior and lateral pelvic stability

Reflex[™] Foam



- Perforated in sacral-ischial area
 - Creates lower stiffness under ITs
 - Simulates grey foam of the Evolution[™]

Air-Release Device

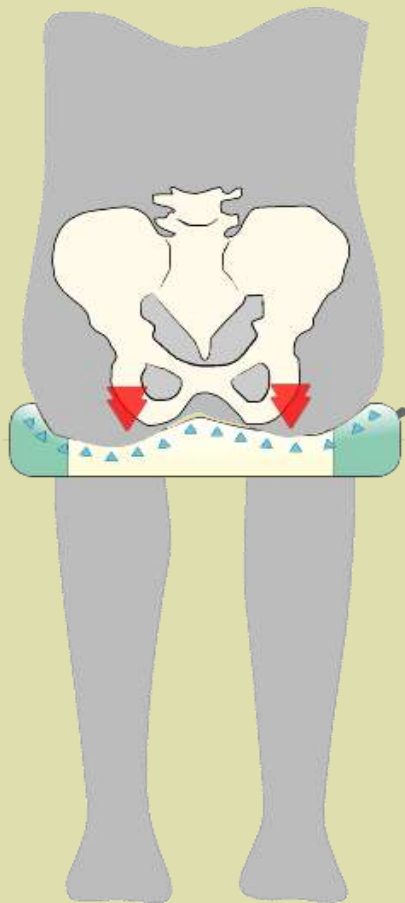

- Built in air-release device
- Automatically releases a portion of air when client sits on cushion
- Allows a fixed amount of immersion
- Prevents bottoming out
- Automatically re-inflates when weight is taken off cushion

Fabric



- Bonded to foam
- Puncture resistant
- Water resistant
- Provides strength and durability
- Grey polyester

Cover



- Tailored
- Removable
- Incontinence resistant
- Non-slip material on bottom
- Hook and loop for added security
- Machine washable
- No reticulated foam
- No zipper, uses Velcro™ for closure

Clients to Consider:

- Ideal for clients that need moderate skin protection and have only mild positioning needs
- Part time ambulators
- A user that wants more comfort than a general use cushion
- Maintenance free feature good for:
 - Understaffed institutions
 - Caregivers
 - Clients with cognitive issues
 - Clients with poor fine motor skills

